

Fume cupboard range Epure

100% French manufacturing
Optimised maintenance and servicing
Optional electric sash
Available in 1200, 1500 and 1800 mm width

**PERMANENT
STOCK**

- Using **6mm laminated glass for the front panel**, it provides complete visualization of the internal volume.
- Worktop with airflow profile forming an internal retention.
- Structure made of 19 mm class B-s1 d0 melamine board, resistant to humid conditions.
Side panels equipped with a thermoformed Acrovyn profile in a **customizable color** (optional) on which the control box and the power sockets can be fixed.
- Fume cupboards sash composed of an **anodized aluminum frame** and a 6 mm laminated glass (or polycarbonate) can be opened by a system of toothed belts with a 15 mm steel core and counterweights positioned on the front side.
- Handle in anodized aluminum positioned at the bottom to optimize air flow. **Stop system at a height of 400 mm or 500 mm** with manual unlocking (Motorised, optional).
- Back plenum, high-density polyethylene support studs with fixings for **Labomeca® mounting**.
- Ceiling in expanded PVC with two explosion-proof vents and a 250 mm diameter hole for air extraction.
- Access to the upper technical box by **clipped front panel** that can be removed without tools.
- **Waterproof LED lighting** positioned outside the suction volume, providing a minimum of 400 Lux lighting.
- Benchtop fume cupboard (freestanding, optional).
- **Flow controlled thanks to an air speed controller** with regulation (visual and audible alarm controlling containment) as well as a frequency converter.
- As an option, the **façade can be equipped with an automatic inactivity controller** for opening and closing, making it possible to secure the fume cupboard by closing its mobile facade after detection of inactivity (**without specific motorization kit**).



NF EN 14175-3
NF X 15-206

Width available	1200, 1500 or 1800 mm
Usable depth	900 mm
Height	1600 mm


Technical specifications fume cupboard Epure 120

Sash open	Air flow	Air speed	Pressure loss
400 mm	460 m ³ /h	0.35 m/s	7.4 Pa
500 mm	577 m ³ /h	0.35 m/s	11.6 Pa

Fume cupboard Epure 120

Supply from ceiling

SPECIFICATION MODEL EPURE 120
SACH OPEN HEIGHT: 400 mm/500 mm

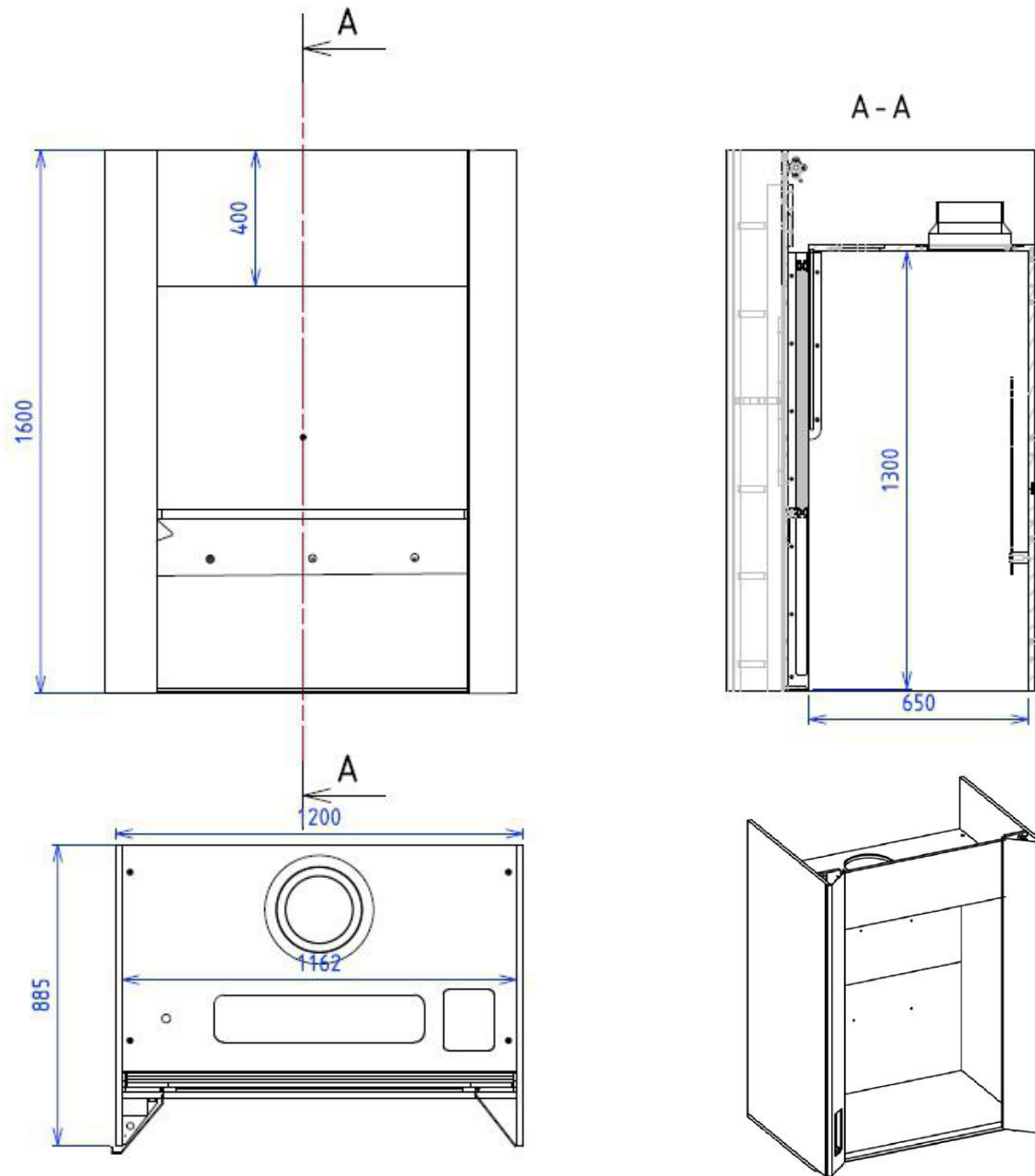

External Dimensions W1200 mm X D 650 mm X H1300 mm
Internal Work Area Dimensions W 1162 mm x D 650 mm x H 1300 mm

EPURE_SORB120	DATE 2022/05/31	PAGE 1/2
 LABOVER	Parc Euromédecine - 20 rue Robert Koch BP 27253 - 34086 Montpellier Cedex 4 Tél. +33 (0)4 67 10 59 00 commercial@labover.fr - www.labover.fr	